

Easy To Use

Simple 2 - step procedure

Perform automatic measurements simply by inserting a tube (for samples), a tip (for sampling and injection), and a cartridge (reagents), and press the START key!

Workflow

1 Insert materials.



③ CARTRIDGE
① TUBE
② TIP

2 Press the START key.



Approx.
10 min

Results

No need to add reagents.

All reagents and materials are already prepared inside the cartridge.

Automated calibration for reagents using attached QR code

All calibrations and measurements are conducted automatically based on the QR code attached to the bottom of the cartridge, which eliminates the need for lot calibration.



"QR code" is a registered trademark of DENSO WAVE INC.

Flexibility in selecting the most appropriate sample tube (1.5 mL or 0.5 mL)

1.5 mL or 0.5 mL sample tube can be used, allowing flexibility when collecting blood samples in accordance with the measurement situation.

Automatic Dilution Function

AU10V has an automatic dilution function for v-COR measurements to accommodate serum samples with high cortisol levels. (Plasma samples are to be measured in dilution mode regardless of cortisol level.)


Analyzer Specifications

Dimensions (W x D x H)	238 x 326 x 310 [mm]
Weight	Approximately 7 kg
Operating temperature	15°C - 30°C (59°F - 86°F)
Setup time	3 min
Required sample volume	basic test : 100 µL dilution test : 20 µL
Data memory	100 tests
Cartridge information	Information about test parameters is embedded in the QR code attached to the cartridge.

Reagent Cartridges Specifications

Types of sample	v-T4 : Canine and Feline Plasma/Serum vc-TSH : Canine Plasma/Serum v-COR : Canine Plasma (dilution) /Serum v-BA : Canine and Feline Plasma/Serum v-PRG : Canine Plasma/Serum
Dynamic Range	v-T4 : 0.50 - 8.00 µg/dL (6.4 - 103.0 nmol/L) vc-TSH : 0.25 - 5.00 ng/mL v-COR : 1.0 - 30.0 µg/dL (27.6 - 828.0 nmol/L) 1.0 - 50.0 µg/dL (27.6 - 1380.0 nmol/L) (with dilution) v-BA : 2.0 - 150.0 µmol/L v-PRG : 0.20 - 40.00 ng/mL (0.64 - 127.20 nmol/L)
Measurement time	Approximately 10 min (without dilution)
Package Unit	10 cartridges/box
Reagent Storage	2°C - 8°C (35.6°F - 46.4°F)

System Configuration Diagram

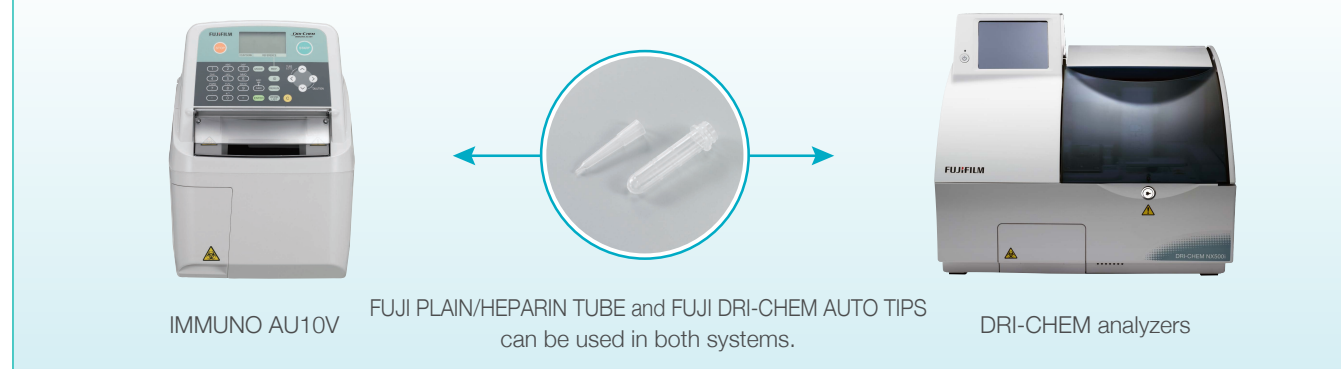

Options

Sample barcode reader




Printer




DIMENSIONS


Top

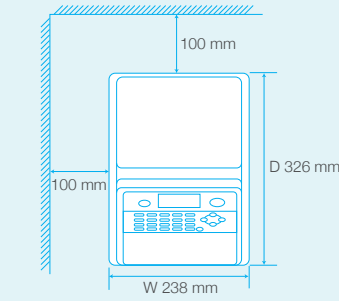

Front

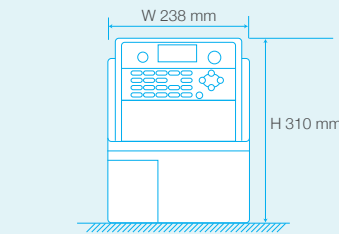

Side

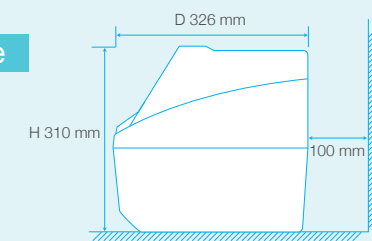

FUJIFILM
Value from Innovation

DRI-CHEM IMMUNO AU10V

Automated Fluorescence Immunoassay Analyzer



NEW Progesterone Test

Innovation in Immunodiagnostic Testing

FUJIFILM

FUJIFILM Corporation

26-30, NISHIAZABU 2-CHOME, MINATO-KU, TOKYO 106-8620, JAPAN
<http://www.fujifilm.com/products/medical/>

CE

Ref. No. FDCAU10V R2 (18-06-F1079-F8598) Printed in Japan ©2011 FUJIFILM Corporation

DRI-CHEM IMMUNO AU10V (Product: FUJI DRI-CHEM IMMUNO AU10V)

This equipment complies with FDA performance standards for laser products except for deviations pursuant to Laser Notice No. 50, dated (June 24, 2007). This equipment is a Class 1 laser product (designed according to IEC60825-1:2014 / EN60825-1:2014)

Key Features

FUJIFILM's new Surface Plasmon enhanced Fluorescence (SPF) technology has brought about a compact sized Fluorometry Immunodiagnostic Analyzer of high reliability and accuracy. Ideal for in-clinic immunodiagnostic testing.


New test parameter: Progesterone (v-PRG)

- Easy to use
- Rapid measurement (approximately 10min)
- Serum and plasma samples are measurable in all tests
- Feline samples are measurable in v-T4 and v-BA
- The only bench top analyzer with TSH test

Reagent Cartridges Specifications

Test Parameter	Application
NEW v-PRG (Progesterone)	Reproduction management
v-T4 (Thyroxine (T4))	Hypo/Hyperthyroidism diagnosis
vc-TSH (Thyroid-Stimulating Hormone)	Hypothyroidism diagnosis
v-COR (Cortisol)	Adisson's Disease or Cushing's syndrome diagnosis
v-BA (Bile Acids)	Liver-related diseases diagnosis

Real time diagnosis and monitoring for hormones are possible now.

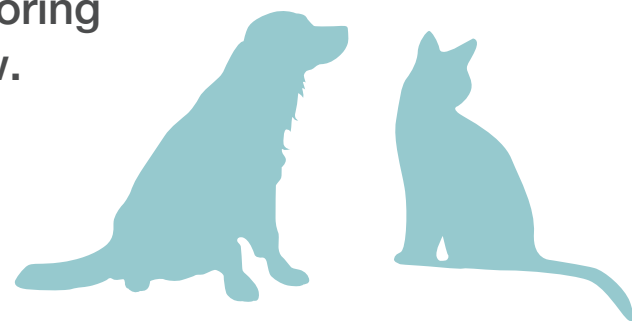

DRI-CHEM IMMUNO AU10V

A new in-clinic Immunodiagnostic Analyzer realized by FUJIFILM's cutting edge technologies




FUJIFILM DRI-CHEM IMMUNO AU CARTRIDGE v-PRG



FUJIFILM DRI-CHEM IMMUNO AU CARTRIDGE v-T4



FUJIFILM DRI-CHEM IMMUNO AU CARTRIDGE vc-TSH

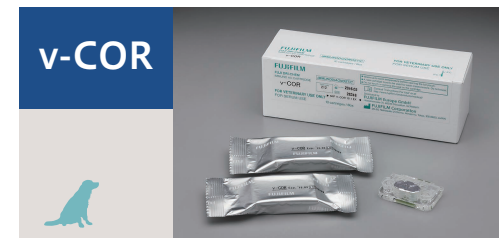

FUJIFILM DRI-CHEM IMMUNO AU CARTRIDGE v-COR





FUJIFILM DRI-CHEM IMMUNO AU CARTRIDGE v-BA



Technology

Principle of the system

- The equipment performs fluorometry using Surface Plasmon enhanced Fluorescence (SPF), a technology based on Surface Plasmon Resonance (SPR), within a small cartridge.
- With the SPF technology, washing for surplus fluorescent beads is not necessary, making the measurement time shorter and the machine more compact. Also since the fluorescence is enhanced by SPR, a low power beam is sufficient enough to generate substantial fluorescence intensity for detection.

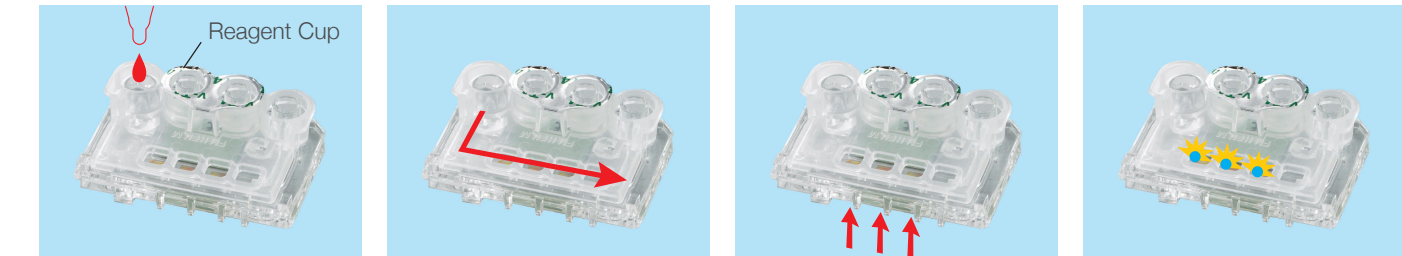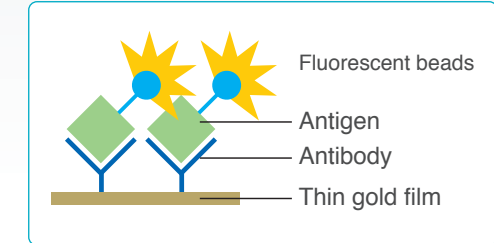

① Sample Injection

The sample, which is mixed with the reagent, will be dispensed into the sample flow channel of the cartridge.

② Binding Reaction

The analyte binds with the antibody, which will be attached to the thin gold film.

③ Laser Beam Irradiation

A SPR is generated upon irradiation by the laser beam.

④ Fluorescence Measurement

The fluorescent intensity, which is commensurate with the level of analyte concentration, is detected.

High Reliability

Excellent correlation with standard laboratory analyzer

AU10V provides reliable and accurate results. It has excellent correlations with its comparison methods in serum samples.

